



VISUAL LANGUAGE FOR CONTACTLESS CARE OPERATIONS

Calm enough to scan. Precise enough to trust.




PRODUCT RULE

Adaptive Care separates resident attention, device health, confidence, freshness, and workflow so staff can decide what to do next.

The brand can speak boldly here. Product UI typography stays compact and restrained.

Product mode	Clinic operations
Visual direction	Clear Signal V3.0
Product type	Restrained Geist Sans
Access floor	WCAG 2.2 AA



01

Foundation

A bright, approachable reference for live care operations.

01

Separate the facts

Risk, device health, confidence, freshness, and workflow are never collapsed into one status.

02

Show uncertainty

Missing or limited evidence is labelled clearly. Last-known values never pretend to be current.

03

Lead with action

One decision area gets one primary next step. Supporting detail stays useful and quiet.

“Which resident needs my attention right now?”

The question every clinic screen should answer within two seconds.

Color

Cobalt guides action. Supporting color stays explicit and useful.

FOUNDATION TOKENS

Canvas

Page background

App shell and section canvas

#FBFAF8

White

Working surface

Cards, forms, and tables

#FFFFFF

Mist 25

Lifted surface

Hover and quiet modules

#FBFCFF

Mist 50

Cool support

Informational section canvas

#F4F7FC

Mist 200

Light border

Card and table boundaries

#DDE4EF

Ink 600

Secondary text

Helper copy and metadata

#536174

Ink 950

Primary text

Headings and decisions

#111827

BRAND AND SUPPORTING TOKENS

Blue 50

Selection wash

Selected rows and active navigation

#EEF4FF

Blue 100

Strong selection

Focused information groups

#DCE9FF

Blue 600

Brand and interaction

Primary actions and focus

#175CD3

Blue 800

Pressed interaction

Pressed primary actions

#123E9A

Sky 400

Informational accent

Device and data information

#55ACFF

Violet 500

Brand accent only

Charts and brand moments, never severity

#7357D8

Mint 400

Positive accent

Positive highlights, not device health

#76D6B1

OPERATIONAL STATUS TOKENS**Positive green**

Operationally positive

Current monitoring and resolved workflow

#147D5A

Watch amber

Review support

Limited coverage or investigation

#A15C00

Risk red

Resident attention

Critical or high attention priority

#C53B30

Unavailable gray

Evidence limitation

Missing, stale, or unavailable evidence

#526172

Color rule Blue means brand or interaction. Sky means information. Violet is a brand accent and never severity. Green, amber, red, and gray carry operational status.

Typography

Compact, direct language built for repeated scanning.

Display · 32 / 38 · 650

Care, without the noise.

Page title · 28 / 34 · 650

Residents needing attention

Major heading · 22 / 28 · 650

Morning care review

Section heading · 18 / 24 · 650

Monitoring coverage

Record title · 15 / 21 · 600

Resident B · Room 214

Body · 14 / 20 · 400

Monitoring is active. Latest room evidence arrived 38 seconds ago.

Body strong · 14 / 20 · 550

Review the current evidence before resolving.

Label · 13 / 18 · 550

Workflow state

Metadata · 12 / 17 · 450

Updated 08:42 · Synthetic record

Overline · 11 / 16 · 650

RESIDENT ATTENTION

Mono reading · 12 / 17 · 500

AC-R214-B 08:42:18 +00:38

Space & shape

A four-pixel grid, open spacing, and soft geometry keep information approachable.

04

4px

08

8px

12

12px

16

16px

24

24px

32

32px


48

48px


64

64px

CORNER RADIUS

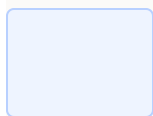**Small**

4px

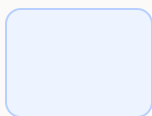**Control**

8px

**Specimen**

20px

**Feature**

24px

DEPTH

Base surface

Hairline border

Popover

Overlap shadow

24px**Default card padding**

Use 16px on narrow screens.

16px**Related component gap**

Use 24px between major modules.

32px**Desktop page gutter**

Reduce to 16px below 768px.

96px**Section spacing**

Keep reference groups distinct.

Actions

Controls are restrained, stable, and explicit about outcome.

BUTTON HIERARCHY

Review resident

Acknowledge

View history

Remove device

...

LOADING AND SUCCESS SPECIMENS

Saving change

✓ Saved

Resolve event

PRIMARY ACTION STATES

RESTING

Review
resident

HOVER VISUAL
SPECIMEN

Review
resident

FOCUS VISUAL
SPECIMEN

Review
resident

DISABLED

Review
resident

ICON ACTION STATES

RESTING

...

HOVER VISUAL
SPECIMEN

...

FOCUS VISUAL
SPECIMEN

...

DISABLED

...

Action rule Disabled actions keep their label and explain the reason nearby. Loading never changes a button's width.

Forms

Visible labels, clear requirements, and local recovery.

Search residents

Search by resident label or room.

Example: Room 214

Search focus state

Visible focus uses the blue interaction ring.

Room 214

Workflow state

Acknowledged



Follow-up timing

☒ This round

☐ Next round

Care note Required

Record what staff observed.

Resident checked.



Record staff observation pending review

Save observation

Add what staff observed before saving.

Status

Six independent axes describe what is known, what is limited, and what staff should do next.

Status rule Attention, monitoring, confidence, freshness, device, and workflow are separate facts. Their colors support the written meaning; they never replace it.

Attention

Resident-facing priority and caregiver urgency.

- ATTENTION**
Critical attention priority
 Highest priority for caregiver review; review the evidence before acting.
- ATTENTION**
High attention priority
 High priority for caregiver review; review the evidence before acting.
- ATTENTION**
Watch attention priority
 Review when practical; this is an attention priority, not a diagnosis.
- ATTENTION**
No attention priority
 No attention priority is indicated.

Monitoring

Whether resident attribution is usable.

- MONITORING**
Monitoring active
 Resident attribution is usable.
- MONITORING**
Monitoring away
 Resident is away; resident-specific baseline learning is paused.
- MONITORING**
Monitoring possible multi-person
 Resident-specific attribution is unavailable; do not guess which person caused a signal.
- MONITORING**
Monitoring paused
 Resident-specific output is unavailable until monitoring resumes.
- MONITORING**
Monitoring calibrating
 Resident-specific output may be unavailable while setup settles.
- MONITORING**
Monitoring unavailable
 No current resident-specific evidence can be shown.

Confidence

Evidence quality and attribution certainty.

- CONFIDENCE
High confidence
Evidence quality supports a high-confidence output.

- CONFIDENCE
Medium confidence
Evidence quality is mixed; review the surrounding context.

- CONFIDENCE
Low confidence
Evidence quality is limited; treat resident-specific output cautiously.

- CONFIDENCE
Confidence unavailable
Evidence quality is unavailable; do not draw a resident-specific conclusion.

Freshness

When evidence was last current; last-known values are not live.

- FRESHNESS
Freshness current
Evidence is current.

- FRESHNESS
Freshness delayed
The latest evidence is delayed; review before treating it as current.

- FRESHNESS
Stale
Last current update: 08:42:18. Current evidence may have changed since the last current update.

- FRESHNESS
Freshness unknown
The last current update is unknown; do not treat values as live.

Device

Room-unit and source health, not resident status.

- DEVICE
Healthy device
Room sources are reporting.

- DEVICE
Degraded device
One or more room sources may be limited.

- DEVICE
Offline device
Room unit is not currently reporting.

- DEVICE
Device maintenance
Room monitoring may be unavailable during device work.

Workflow

Work lifecycle; resolved history remains immutable.

- WORKFLOW
Workflow new
New work item; review event.

- WORKFLOW
Workflow acknowledged
Acknowledged; follow up with the available evidence.

- WORKFLOW
Workflow investigating
Keep the event open while checking.

- WORKFLOW
Workflow resolved
Resolved; this history remains immutable.

Independent facts can coexist

- ATTENTION

High attention priority

High priority for caregiver review; review the evidence before acting.

- MONITORING

Monitoring active

Resident attribution is usable.

- CONFIDENCE

Low confidence

Evidence quality is limited; treat resident-specific output cautiously.

- FRESHNESS

Stale

Last current update: 08:42:18. Current evidence may have changed since the last current update.

- DEVICE

Offline device

Room unit is not currently reporting.

- WORKFLOW

Workflow investigating

Keep the event open while checking.

Next action: Review the event and device connection. These statuses are separate facts; do not treat a stale, offline signal as resident-specific without current supporting evidence.

EVIDENCE TRUTH STACK · SYNTHETIC/TEST-ONLY

Each layer has a different source and responsibility.

01 · SENSOR EVIDENCE

Observed input and data quality

Synthetic/test-only example: room telemetry, source contact, and timestamps are observed inputs. They are not a medical conclusion.

02 · DETERMINISTIC WARNING

Rule-based warning

Synthetic/test-only rule output stays visible and cannot be suppressed by AI. Staff review the evidence and policy context.

03 · AI INTERPRETATION

Possible interpretation

Plain-language interpretation of structured evidence only; it is never a diagnosis and cannot invent facts.

04 · STAFF OBSERVATION

Human-entered context

Record what staff observed or checked. This context supports workflow and does not rewrite immutable history.

Data

Repeated records use tables. Missing values say what they mean.

RESIDENT RECORDS · SYNTHETIC/TEST-ONLY

Synthetic resident monitoring records			
RESIDENT	ATTENTION PRIORITY	EVIDENCE	OPERATIONS
Avery Chen Room 214 ATTENTION REASON No attention priority; monitoring is current	ATTENTION No attention priority	CONFIDENCE High confidence FRESHNESS Freshness current	MONITORING Monitoring active DEVICE Healthy device WORKFLOW Workflow acknowledged Review record
Jordan Lee Room 108 ATTENTION REASON No attention priority; a current reading is available HOVER EXAMPLE	ATTENTION No attention priority	CONFIDENCE High confidence FRESHNESS Freshness current	MONITORING Monitoring active DEVICE Healthy device WORKFLOW Workflow acknowledged Open record
Sam Rivera Room 302 ATTENTION REASON Multiple people may be present in the room SELECTED RECORD	ATTENTION No attention priority	CONFIDENCE Confidence unavailable FRESHNESS Freshness unknown	MONITORING Monitoring possible multi-person Resident attribution unavailable; do not guess which resident caused this signal. DEVICE Healthy device WORKFLOW Workflow new Review record

RESIDENT	ATTENTION PRIORITY	EVIDENCE	OPERATIONS
Priya Shah Room 220 ATTENTION REASON Resident-away period needs a coverage check	ATTENTION Watch attention priority	CONFIDENCE Confidence unavailable FRESHNESS Stale last current update 08:38:12	MONITORING Monitoring away DEVICE Healthy device WORKFLOW Workflow acknowledged Review record
Mateo Brooks Room 305 ATTENTION REASON Unexpected movement needs review	ATTENTION Critical attention priority	CONFIDENCE Low confidence FRESHNESS Freshness delayed	MONITORING Monitoring active DEVICE Degraded device WORKFLOW Workflow investigating Review record

Responsive rule Desktop uses a scannable table; narrow screens use ordered records so the resident, reason, evidence, workflow, and action remain visible without horizontal panning.

Care patterns

Operational components keep priority, truth, workflow, and action together.

ATTENTION QUEUE ITEM

Resident B Room 214

12 min open

ATTENTION REASON

Unexpected movement needs review

- ATTENTION

High attention priority

High priority for caregiver review; review the evidence before acting.

- MONITORING

Monitoring active

Resident attribution is usable.

- CONFIDENCE

Low confidence

Evidence quality is limited; treat resident-specific output cautiously.

- FRESHNESS

Freshness current

Evidence is current.

- DEVICE

Healthy device

Room sources are reporting.

- WORKFLOW

Workflow investigating

Keep the event open while checking.

OBSERVED CONTEXT

Synthetic/test-only example. Review the available evidence before deciding what to do next.

[Review resident](#)

► [Device details for Resident B](#)

RESIDENT STATUS

Monitoring current

High confidence. Latest evidence arrived 38 seconds ago.

Keep risk, confidence, and freshness separate.

DEVICE STATUS

Healthy device

Radar, thermal, and Wi-Fi sources are reporting.

Device health is separate from resident attention and confidence.

ATTENTION HIERARCHY

- 1 High attention priority
- 2 Review workflow state
- 3 Device and data limits

This synthetic example does not identify a medical cause.

OPERATIONAL STATES · SYNTHETIC/TEST-ONLY

One reusable catalog keeps known facts, limitations, why they matter, and the allowed next action together. Color supports the words; it never carries meaning alone.

MONITORING

Info

Calibrating

Provenance: Setup and baseline state

What is known

The room setup is settling and calibration work is in progress.

What is unavailable or limited

Resident-specific output may be limited; an established baseline is not ready yet.

Why it matters

Early readings should not be treated as a stable picture of this resident's routine.

Allowed next action

Wait for calibration to complete, then review the resulting baseline status.

MONITORING

Caution

Partial baseline

Provenance: Resident baseline state

What is known

Some resident-specific baseline dimensions have been established.

What is unavailable or limited

Other dimensions do not have enough eligible history for a complete comparison.

Why it matters

An anomaly comparison may be less informative while the baseline is partial.

Allowed next action

Use the available evidence cautiously and allow eligible monitoring history to accumulate.

MONITORING

Limited

Resident away

Provenance: Occupancy and monitoring state

What is known

The resident-away state is active for the assigned room.

What is unavailable or limited

Resident-specific monitoring is paused while the resident is away.

Why it matters

Resident-specific baseline learning is paused so absence is not learned as normal behavior.

Allowed next action

Confirm the coverage context when needed and wait for a return transition before resident-specific review.

MONITORING

Info

Resident returned / return transition

Provenance: Occupancy transition

What is known

The room has a new return transition after an away period.

What is unavailable or limited

A return does not by itself prove that the resident is safe or that all resident-specific evidence is current.

Why it matters

Returning is an awareness update, not a clinical conclusion; learning resumes only when eligible.

Allowed next action

Review current evidence and let the monitoring state settle before making a resident-specific decision.

MONITORING

Limited

Possible multi-person presence

Provenance: Occupancy ambiguity state

What is known

Room telemetry suggests that more than one person may be present.

What is unavailable or limited

Resident-specific attribution is unavailable; do not guess which person caused a signal.

Why it matters

Resident-specific baseline learning is paused while attribution is ambiguous.

Allowed next action

Confirm occupancy context before using resident-specific output or recording a resident-specific conclusion.

DEVICE

Caution

Device degraded

Provenance: Device-health state

What is known

The room device is reporting, but one or more sources may be limited.

What is unavailable or limited

Some current sensor evidence may be missing or lower quality.

Why it matters

Device health is separate from resident attention; limited sources can reduce confidence without creating a resident conclusion.

Allowed next action

Review device and source health, then use only the evidence that is currently available.

DEVICE

Limited

Device offline

Provenance: Device-health state

What is known

The room unit is not currently reporting to the monitoring service.

What is unavailable or limited

Current room and resident-specific evidence is unavailable until reporting resumes.

Why it matters

Device health is separate from resident attention; an offline unit is not evidence that the resident is safe or unsafe.

Allowed next action

Check the room unit's connection and follow the device recovery procedure.

TRANSPORT

Caution

Buffering

Provenance: Edge transport state

What is known

The edge unit is holding compact telemetry while upload is unavailable or delayed.

What is unavailable or limited

The held telemetry is not current evidence; the interface does not fabricate current data.

Why it matters

A growing transport delay can make resident-specific output stale even when the device is powered.

Allowed next action

Review connection health and wait for confirmed ingestion before treating a reading as current.

TRANSPORT

Caution

Retrying

Provenance: Transport recovery state

What is known

A failed upload or operation is being attempted again.

What is unavailable or limited

Current data may be delayed; retrying does not fabricate current data or confirm success.

Why it matters

The next decision should use confirmed evidence, not an assumed successful retry.

Allowed next action

Allow the retry to finish, then verify the last confirmed update before acting.

ASSIGNMENT

Limited

Missing room/resident assignment

Provenance: Room and resident configuration

What is known

The device or room does not have a complete authorized resident assignment.

What is unavailable or limited

Resident-specific attribution and resident baseline comparisons are unavailable.

Why it matters

Room telemetry must not be presented as belonging to a resident without an authorized assignment.

Allowed next action

Complete or verify the room and resident assignment before using resident-specific output.

ANOMALY

Caution

Unknown anomaly

Provenance: Deterministic anomaly output

What is known

A meaningful deviation was detected, but it does not match a supported event family.

What is unavailable or limited

The system cannot name a cause from this evidence and makes no diagnosis.

Why it matters

Keeping an unknown anomaly general avoids forcing uncertain evidence into a misleading label.

Allowed next action

Review the structured evidence and observed context; record what staff can verify.

INTERPRETATION

Info

AI interpretation pending

Provenance: AI interpretation stage

What is known

A structured evidence packet is awaiting an optional plain-language interpretation.

What is unavailable or limited

The interpretation is not ready; the deterministic warning remains visible and actionable.

Why it matters

AI adds context after evidence exists; it does not monitor telemetry or decide whether a warning is real.

Allowed next action

Review the deterministic evidence now and wait for interpretation only if additional context is useful.

INTERPRETATION

Limited

AI interpretation unavailable

Provenance: AI interpretation stage

What is known

The deterministic evidence and warning decision remain available.

What is unavailable or limited

The optional AI interpretation is unavailable; it is never a diagnosis.

Why it matters

The deterministic warning remains visible and actionable; a missing interpretation cannot suppress it or remove the work item.

Allowed next action

Review the deterministic warning and structured evidence; continue the workflow without AI interpretation.

FEEDBACK

Info

Loading

Provenance: Interface feedback state

What is known

The requested operational content has not finished loading.

What is unavailable or limited

Loading state; current operational details are not available yet.

Why it matters

A placeholder must not look like a current reading or an empty result.

Allowed next action

Wait for the request to complete before making an operational decision.



FEEDBACK

Neutral

Genuine empty

Provenance: Operational collection state

What is known

No items match the requested collection because there are genuinely none to show.

What is unavailable or limited

Nothing is unavailable for this empty collection; monitoring continues where configured.

Why it matters

An explicit empty result prevents staff from mistaking a clear queue for a loading or error state.

Allowed next action

Continue to the next review area or wait for new work to arrive.

FEEDBACK

Neutral

Filtered empty

Provenance: Collection and filter state

What is known

The collection has content, but the active filters match no items.

What is unavailable or limited

Items outside the active filter are not shown in this view.

Why it matters

A filtered empty result is different from a genuinely empty collection.

Allowed next action

Review or clear the filters to see the available collection.

FEEDBACK

Limited

Stale data

Provenance: Evidence freshness state

What is known

A last-known update exists, but it is no longer current.

What is unavailable or limited

The current resident or device condition is unavailable; last-known values are not live.

Why it matters

Stale evidence can change the priority of a decision and must not be shown with fake precision.

Allowed next action

Refresh or check the transport and device state before relying on the last-known value.

FEEDBACK

Risk

Save failure

Provenance: User action result

What is known

The requested change was not confirmed as saved.

What is unavailable or limited

The record may still contain its previous value; the local change is not confirmed.

Why it matters

Unconfirmed workflow changes must not be presented as completed care documentation.

Allowed next action

Retry the save and verify the confirmation before leaving this workflow.

FEEDBACK

Caution

Conflicting update

Provenance: Concurrent update protection

What is known

Another update changed the record before this change could be applied.

What is unavailable or limited

The local change was not applied, and the current record needs to be reloaded.

Why it matters

Conflict recovery prevents one person's newer observation from being silently overwritten.

Allowed next action

Reload the current record, compare the changes, and reapply only what is still correct.

WORKFLOW

Risk

Overdue work

Provenance: Care workflow state

What is known

An open work item passed its expected follow-up time.

What is unavailable or limited

The item remains open until an authorized workflow action resolves it.

Why it matters

High and critical events never silently expire; overdue status needs explicit follow-up.

Allowed next action

Review the available evidence and continue, escalate, or resolve through the documented workflow.

WORKFLOW

Positive

Resolved read-only

Provenance: Immutable event history

What is known

The work item has a recorded resolution and is shown as history.

What is unavailable or limited

Resolved history cannot be edited or reopened in place.

Why it matters

Resolved history remains immutable so the record preserves what happened and when.

Allowed next action

Review the recorded outcome; create a new linked event if the same concern occurs again.

WORKFLOW

Caution

Recurrence / new linked event

Provenance: Event relationship and workflow state

What is known

A new event is linked to an earlier resolved event as a recurrence.

What is unavailable or limited

The earlier resolved event remains read-only and is not reopened.

Why it matters

A new linked event preserves separate timelines, ownership, and evidence for each occurrence.

Allowed next action

Review the new event as current work while using the linked history for context.

WORKFLOW SPECIMEN · SYNTHETIC/TEST-ONLY

Event lifecycle

Open work progresses explicitly. Resolved history is immutable; a recurrence is a new linked event.

01

Open · new

A new work item is ready for evidence review.

Review the available evidence before acknowledging or investigating.

02

Acknowledged

The work item is seen and remains open for follow-up.

Continue with the documented check; acknowledgement does not resolve the item.

03

Investigating

Staff are checking the evidence and observed context.

Keep the item open until the check has a recorded outcome.

04

Overdue

Expected follow-up time passed while the item remained open.

Take explicit follow-up action; high and critical events never silently expire.

05

Resolved read-only

The event has a recorded outcome and remains immutable history.

Read the outcome; never reopen resolved history in place.

06

Recurrence / new linked event

A new event links to the earlier resolved event for context.

Work the new event as current work while preserving the prior event unchanged.

Watch items may auto-close into history only when that behavior is documented; high and critical events never silently expire.

INTERFACE FEEDBACK · SYNTHETIC/TEST-ONLY

System feedback

Skeleton loading, genuine empty, filtered empty, stale refresh, retryable save failure, and conflict recovery explain whether content is loading, absent, old, or needs recovery.

FEEDBACK

Info

Loading

Provenance: Interface feedback state

What is known

The requested operational content has not finished loading.

What is unavailable or limited


Loading state; current operational details are not available yet.

Why it matters

A placeholder must not look like a current reading or an empty result.

Allowed next action

Wait for the request to complete before making an operational decision.

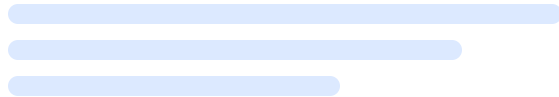

FEEDBACK

Neutral

Genuine empty

Provenance: Operational collection state

What is known

No items match the requested collection because there are genuinely none to show.

What is unavailable or limited

Nothing is unavailable for this empty collection; monitoring continues where configured.

Why it matters

An explicit empty result prevents staff from mistaking a clear queue for a loading or error state.

Allowed next action

Continue to the next review area or wait for new work to arrive.

FEEDBACK

Neutral

Filtered empty

Provenance: Collection and filter state

What is known

The collection has content, but the active filters match no items.

What is unavailable or limited

Items outside the active filter are not shown in this view.

Why it matters

A filtered empty result is different from a genuinely empty collection.

Allowed next action

Review or clear the filters to see the available collection.

FEEDBACK

Limited

Stale data

Provenance: Evidence freshness state

What is known

A last-known update exists, but it is no longer current.

What is unavailable or limited

The current resident or device condition is unavailable; last-known values are not live.

Why it matters

Stale evidence can change the priority of a decision and must not be shown with fake precision.

Allowed next action

Refresh or check the transport and device state before relying on the last-known value.

FEEDBACK

Risk

Save failure

Provenance: User action result

What is known

The requested change was not confirmed as saved.

What is unavailable or limited

The record may still contain its previous value; the local change is not confirmed.

Why it matters

Unconfirmed workflow changes must not be presented as completed care documentation.

Allowed next action

Retry the save and verify the confirmation before leaving this workflow.

FEEDBACK

Caution

Conflicting update

Provenance: Concurrent update protection

What is known

Another update changed the record before this change could be applied.

What is unavailable or limited

The local change was not applied, and the current record needs to be reloaded.

Why it matters

Conflict recovery prevents one person's newer observation from being silently overwritten.

Allowed next action

Reload the current record, compare the changes, and reapply only what is still correct.